

CIVCO Technical Data Sheet

C-Qual M™ Breastboard featuring the Monarch™ Overhead Arm Positioner




Product Name	C-Qual M™ Breastboard	
Product Number	MTM400	C-Qual M Supine Breastboard with Monarch Overhead Arm Positioner
	OAP100	Monarch Overhead Arm Positioner for use with C-Qual™
	OAP110	Monarch Overhead Arm Positioner for use with Lok-Bar™
General Description	The C-Qual M assists in providing comfortable, reproducible patient positioning for supine breast treatment. The board features an enhanced bottom stop. The board integrates with the Monarch Overhead Arm Positioner to provide an efficient approach for arms-up positioning, allowing for proper superior/inferior adjustment of the board to the patient which helps reduce patient slide.	
Indications for use	Breast cancer patients unable to lie prone. Smaller-breasted patients. Need for clearer visualization of medial field edge. Enhanced CT-gantry clearance.	
MRI Compatible	OAP110 MR Conditional. Product poses no known hazard in 3.0T MR environments.	
Integration	3-Pin Indexable to Couchtop	
FDA Clearance	Yes	
CE Marked (EU DOC)	Yes	
Other Country Registrations	Available upon request	

Technical Details

Primary Materials	Kevlar, Black Acetal, Grey Kydex, Black ABS
Weight (approx)	10.89 kg / 24 lbs
Dimensions (mm)	1333 L x 594 W x 261 H (min.)
Attenuation	1.87% at 6MV thru the middle of the back support

Technical Drawing

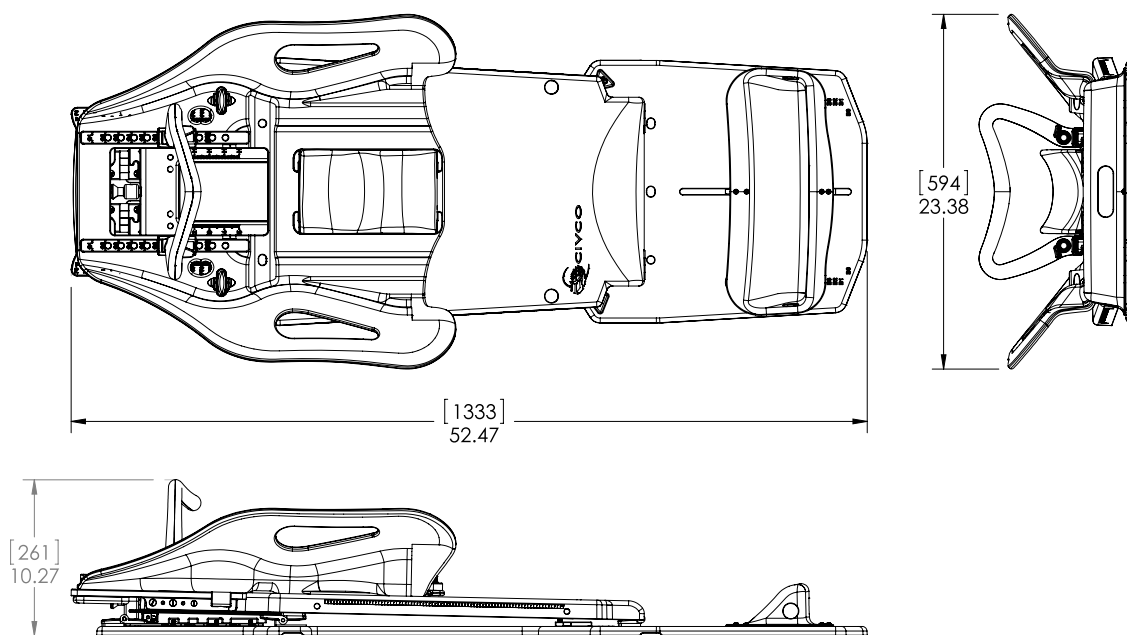

Country of Manufacture	USA
-------------------------------	-----

Packaging

Unit Pack	Each
Units / Box	One device - may have several component parts in box
Shelf Life	n/a
Warranty	See CIVCO website for Product Warranty Statement
Sterilization	Non-Sterile
Launch Date	October 2015

Accessories

System Customization	Compatible with: Posicast® Thermoplastics Posifix® Headrests Breast Support Thermoplastics Vac-Lok™ Cushions C-Qual Storage: 013156
-----------------------------	--

Committed to Positively Impacting Patient Care

Global Offices
Orange City, Iowa, USA Rotterdam, Netherlands info@zivco.com
800.842.8688 | +1.712.737.8688 +31 (0) 10.3036600 www.zivco.com/ro



Radiation Oncology