Body Pro-Lok Rails-Only SBRT Immobilization for use with Exact® & IGRT-style couchtops

The Body Pro-Lok Rails-Only System provides a method for utilizing Body Pro-Lok Bridges, Respiratory Belt and accessories without the full Body Pro-Lok Platform or CIVCO's Universal Couchtop.

Features

- Exact Indexing
- Toggle and clamp design provides secure couchtop attachment
- Reduced attenuation when compared with the use of the Body Pro-Lok Platform
- Rails-Only T-Pin Lok-Bars provide compatibility with indexable products






53cm, Varian Exact	MTSBRT060
Kit includes Exact Rails, Type 1B Bridge, Type 2 Bridge, Respiratory Plate,	MTSBRT060K
Respiratory Belt, 4x T-Pin Lok-Bars™, Kneefix™ 3, Feetfix™ 3, Clam-Lok ™ Cushion	

Secure Couchtop Attachment


The system features four lateral clamps that toggle between an Exact or IGRT profile, allowing a secure fit to each table.

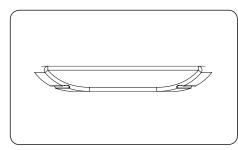

IGRT-style profile, couchtop and lateral clamp - Linac





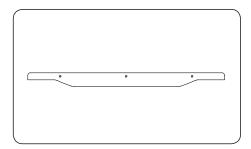

Exact profile, couchtop and lateral clamp - CT

See back for Rails-Only customization







Body Pro-Lok™ Rails-Only Customization

Body Pro-Lok Bridges provide additional immobilization from above.



Type 1B Bridge (front)	MTSBRT039
Type 2 Bridge (middle)	MTSBRT007
Type 3 Bridge (back)	MTSBRT031

The **Clam-Lok™ Cushion** (8") forms to comfortably enhance immobilization from above.



MTSBRT202

The **Feetfix™ 3 Cushion** contributes to better stability and reproducibility in the supine position.



.....

T-Pin Lok-Bars™ allow most indexable positioning aids to be indexed to the Body Pro-Lok Rails-Only System.



Rails-Only Two-Pin Lok-Bar MTSBRT055

Global Sales Office Orange City, Iowa, USA 800.842.8688 | +1 712.737.8688

info@CivcoRT.com www.CivcoRT.com

COPYRIGHT © 2017. CIVCO IS A REGISTERED TRADEMARK OF CIVCO MEDICAL SOLUTIONS. BODY PRO-LOK ONE, ONEBRIDGE, ONEPLATFORM, AND UNIVERSAL COUCHTOP ARE TRADEMARKS OF CIVCO. ALL OTHER TRADEMARKS ARE PROPERTY OF THEIR RESPECTIVE OWNERS. ALL PRODUCTS MAY NOT BE LICENSED IN ACCORDANCE WITH CANADIAN LAW. 2015P1031 REV C

The **Respiratory Plate** places pressure at the level of the diaphragm to assist in restricting respiratory movement.



MTSBRT003

T-Shape Vac-Lok™ Cushion provides positioning with the arms up when used with a Wing Board™, while cradling the torso region.



MTVL0017

The **Extended Wing Board™** is 22.5" (57cm) wide, constructed of durable, lightweight ABS and can be used with T- or U-Grip handles.



with T-Grip Handles with U-Grip Handles



Patient Handle Patient Handle (mirror)

their hands.

The Body Pro-Lok Respiratory Belt

for SBRT treatments.

supine position.

Elevation Blocks

provides pneumatic compression and immobilization to the abdominal region

The Kneefix™ 3 Cushion contributes to

better stability and reproducibility in the

Patient Handles clamp to the platform or

Universal Couchtop[™], providing a stable

place for immobilized patients to place

MTSBRT005 MTSBRT0051

MTSBRT004

301052

301055

