

TotalRT All in One Platform

CIVCO RT Technical Data Sheet




Product Name	TotalRT All in One Platform
Product Number	TRT100
General Description	The TotalRT All in One Platform features Prodigy™ 2 indexing style and profile to quickly adapt CIVCO RT's positioning and immobilization solutions to a patient's setup needs and treatment type. The platform is indexed on a couchtop creating an affordable all in one solution and allowing for seamless transition between patient treatments.
Indications for use	The device is intended to aid in support and positioning of adult and pediatric patients undergoing radiation therapy including electron, photo, and proton treatments. The device is also used during image acquisition to support treatment planning.
Contraindications	None
MRI Compatible	Device is MR Safe
Integration	Connects to couchtop via two or three-pin Lok-Bar™
FDA Clearance	Yes
CE Marked (EU DOC)	Yes
Other Country Registrations	Available upon request
Country of Manufacture	USA


Technical Details

Primary Materials	Glass Fiber, Foam Core			
Weight (approx)	7.3 kg / 16 lbs			
Dimensions (mm)	2013.4 L x 619.71 W x 34.93 H (min.)			
Dosimetric Properties	Location	Field Size (cm)	Attenuation 6 MV photons	Attenuation 18 MV photons
	Middle of platform ●	10 x 10	2.23	1.61
	Head area ●	10 x 10	1.76	1.24
Performance Specifications are for reference only, values may vary based on setup and equipment.				

Technical Drawing

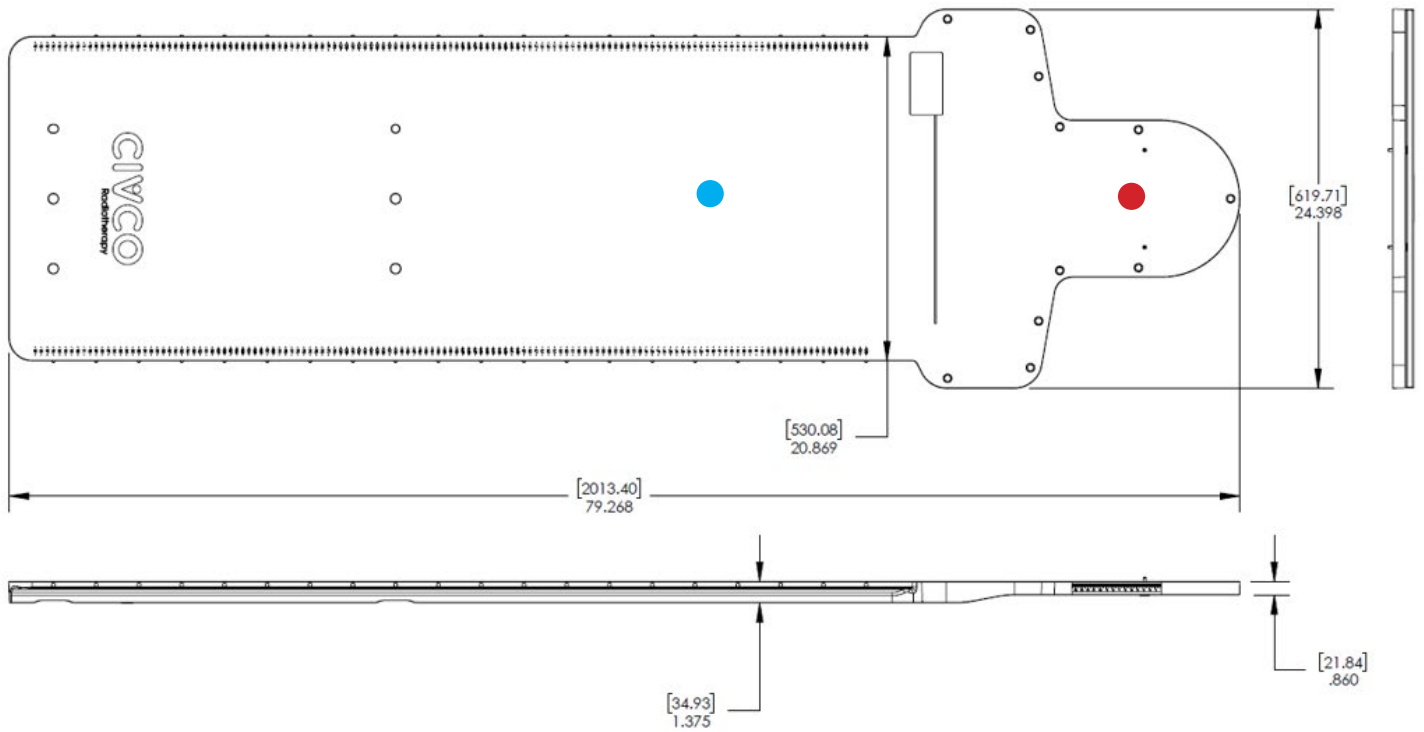

Artifact-Free Treatment Area (Homogeneous Density)

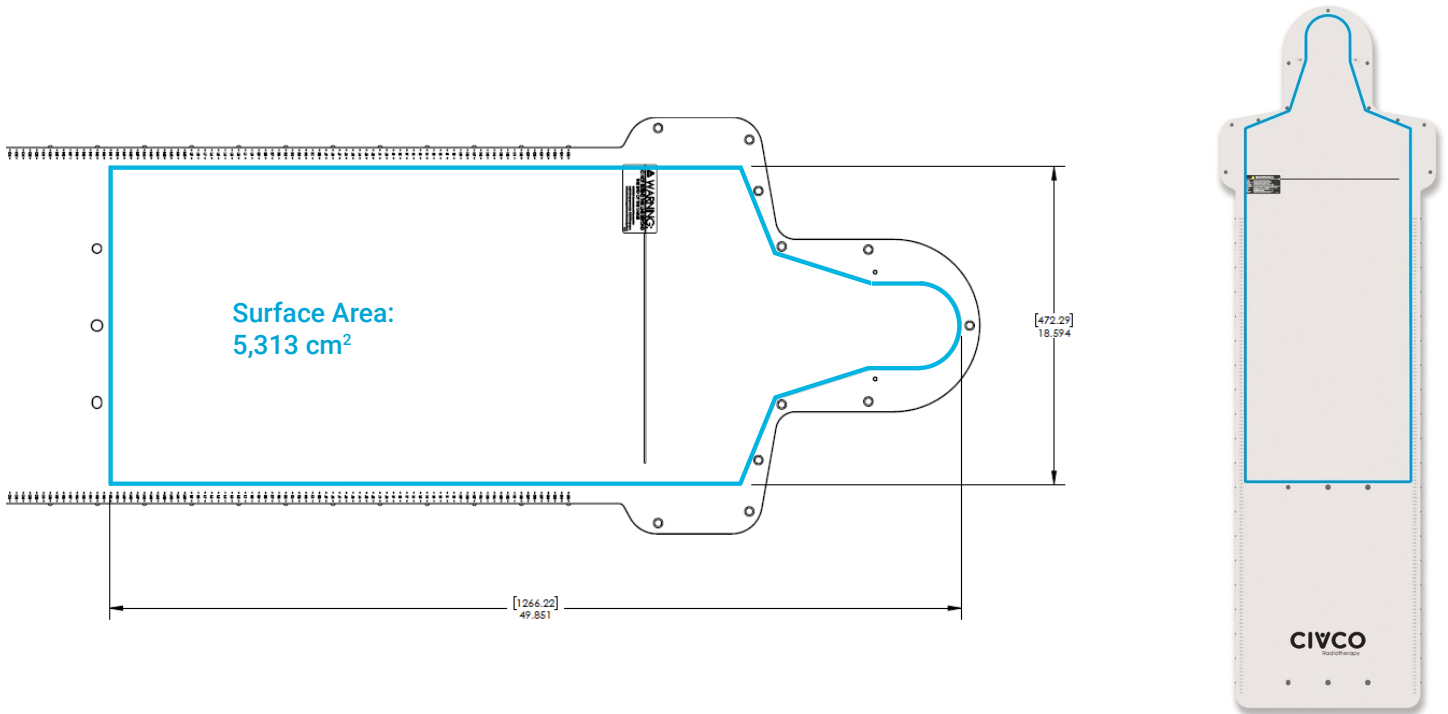

Packaging

Unit Pack	1 Piece
Units / Box	1 Piece
Shelf Life	n/a
Warranty	See CIVCO RT website for Product Warranty Statement
Sterilization	Non-Sterile
Launch Date	October 2021

Accessories

System Customization	CIVCO Radiotherapy positioning and immobilization devices attach directly to platform, or via 2 or 3-Pin Lok-Bars.
-----------------------------	--

Global Sales Office

Orange City, Iowa, USA
800.842.8688 | +1 712.737.8688

info@CivcoRT.com
www.CivcoRT.com

COPYRIGHT © 2022. CIVCO IS A REGISTERED TRADEMARK OF CIVCO MEDICAL SOLUTIONS. CIVCO RADIOTHERAPY, CIVCO RT, LOK-BAR AND PRODIGY ARE TRADEMARKS OF MEDTEC LLC. ALL OTHER TRADEMARKS ARE PROPERTY OF THEIR RESPECTIVE OWNERS. ALL PRODUCTS MAY NOT BE LICENSED IN ACCORDANCE WITH CANADIAN LAW.
2021C1591 REV. B

CIVCO
Radiotherapy™