


Hip & Pelvic Solutions

CIVCO[®]
Radiotherapy
a better way

[Solutions Guide Catalog](#)

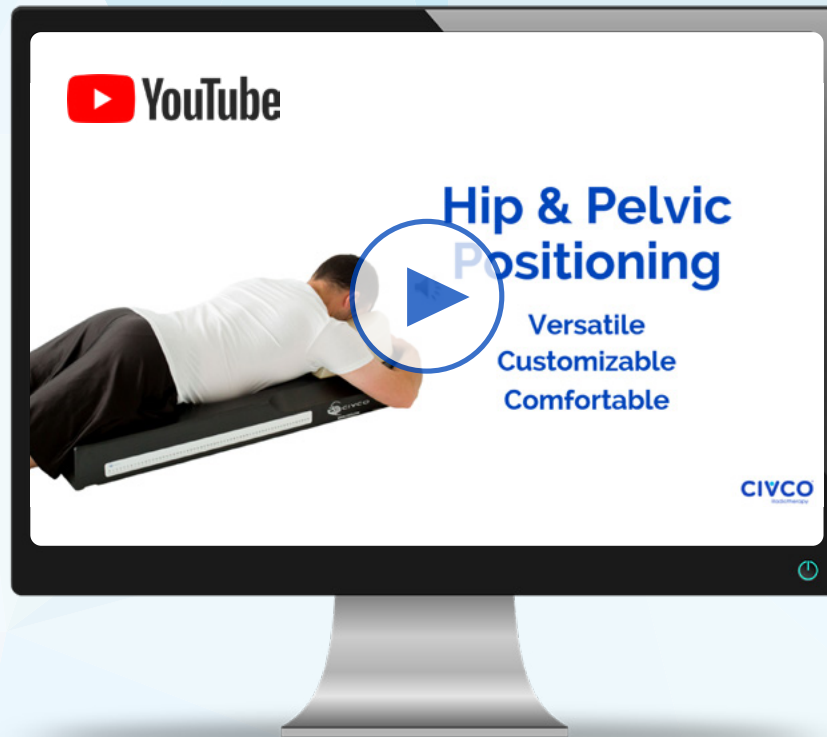

Technical Data Sheets

[Bellyboards](#)

[Combifix™ 3](#)



Sign up for CIVCO Updates!

First Name:

Last Name:

Email:


Organization:

I have read and agree to the [privacy policy](#)



Carbon Fiber & MRSeries™ Bellyboards

The Carbon Fiber Bellyboard assists in reducing the volume of the small bowel that is irradiated in patients undergoing treatment in the pelvic region. The contoured opening allows for displacement of the small bowel, as the central support structure assists in compressing the small bowel forward into the opening.

- Carbon Fiber model features embedded fiducial crosshair to assist as a secondary setup reference
- Lightweight, rigid and indexable
- Lower part of board is adapted to support the patient's upper legs for maximum comfort
- Soft polyethylene padding ensures a comfortable position for the head

125012	Carbon Fiber Bellyboard
126000	MRSeries Bellyboard 
013140	Bellyboard Wall Storage


Caudal Indentation [1] to relieve pressure in the genital area. **Central Space [2]** allows compression to push the small bowel superior toward the **Major Aperture [3]**.

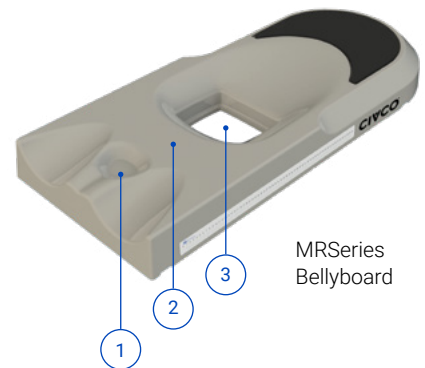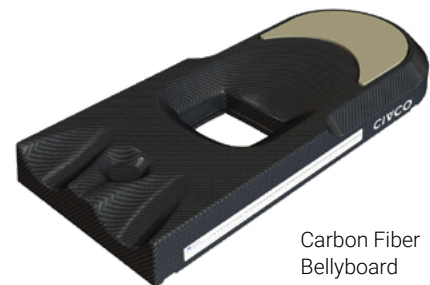

Prone Thorax Support

The Prone Thorax Support is designed to relieve chest compression during treatments in the prone position. It can be used as a stand-alone patient support or as an option on the Bellyboard to allow large patients to be treated in the prone position.

301076	Cream
301102	Black

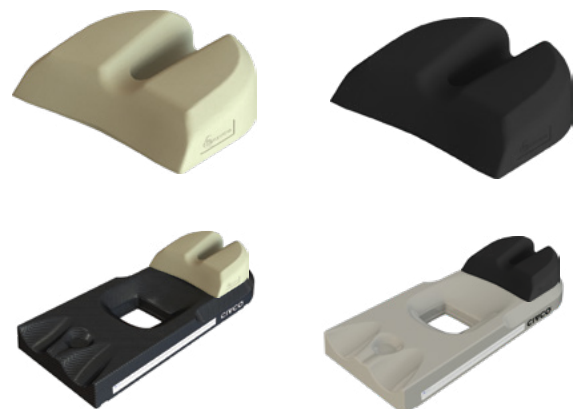

Combifix™ 3 and Feetfix™ Couch Support

The Combifix 3 combines the Kneefix™ and Feetfix, providing enhanced positioning for the pelvic region and other lower extremities.

The Feetfix Couch Support may be used separately and positions the Feetfix at indexed positions, allowing the Feetfix to be tilted at four different angles. Both the Combifix & Feetfix are two and three-pin Lok-Bar indexable.

Shoulder retractors are available to pull the shoulders down during radiotherapy simulation and treatment sessions where opposed lateral fields are used in the neck area. They can attach in five positions on the Feetfix.



Combifix 3 Low



Combifix Storage

151004MR Includes Combifix Baseplate, Kneefix 3, Elevation Blocks (x2), Feetfix 3 & Feetfix Couch Support

151005MR Includes Combifix Baseplate, Low Kneefix 2 Cushion, Elevation Blocks (x2), Feetfix 3 & Feetfix Couch Support

013115 Combifix Wall Storage (not for use in MR environment)

302071MR Feetfix Couch Support Bracket (Cushion not included)

302072 **Shoulder Retractors:** Short (x2) and Long (x2)



Feetfix Couch Support with Cushion



Shoulder Retractors